

LAMP • U



DENTISTRY




Works Brilliantly on Any Devices.


LAMP.U works really well on a wide range of digital devices — from Android or iOS smartphones to compact digital cameras.

Perfect Images Everytime.

LAMP.U is designed with the beginner dental photographer in mind. Ergonomic to use with its chin guard setup, LAMP.U takes perfect standardized picture - everytime.

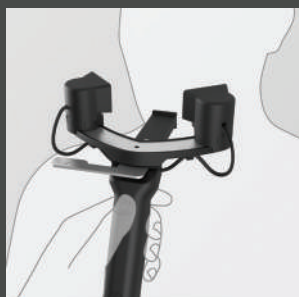


Dental Photography Made Easier Than Ever.

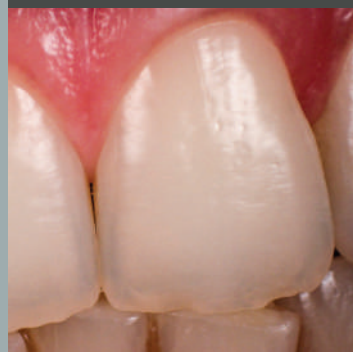

Empowering You.

LAMP.U lets you take standardized dental photos using smartphones and compact cameras that you already own. Now anyone can take professional photos with clarity, crisp details, and realistic color.

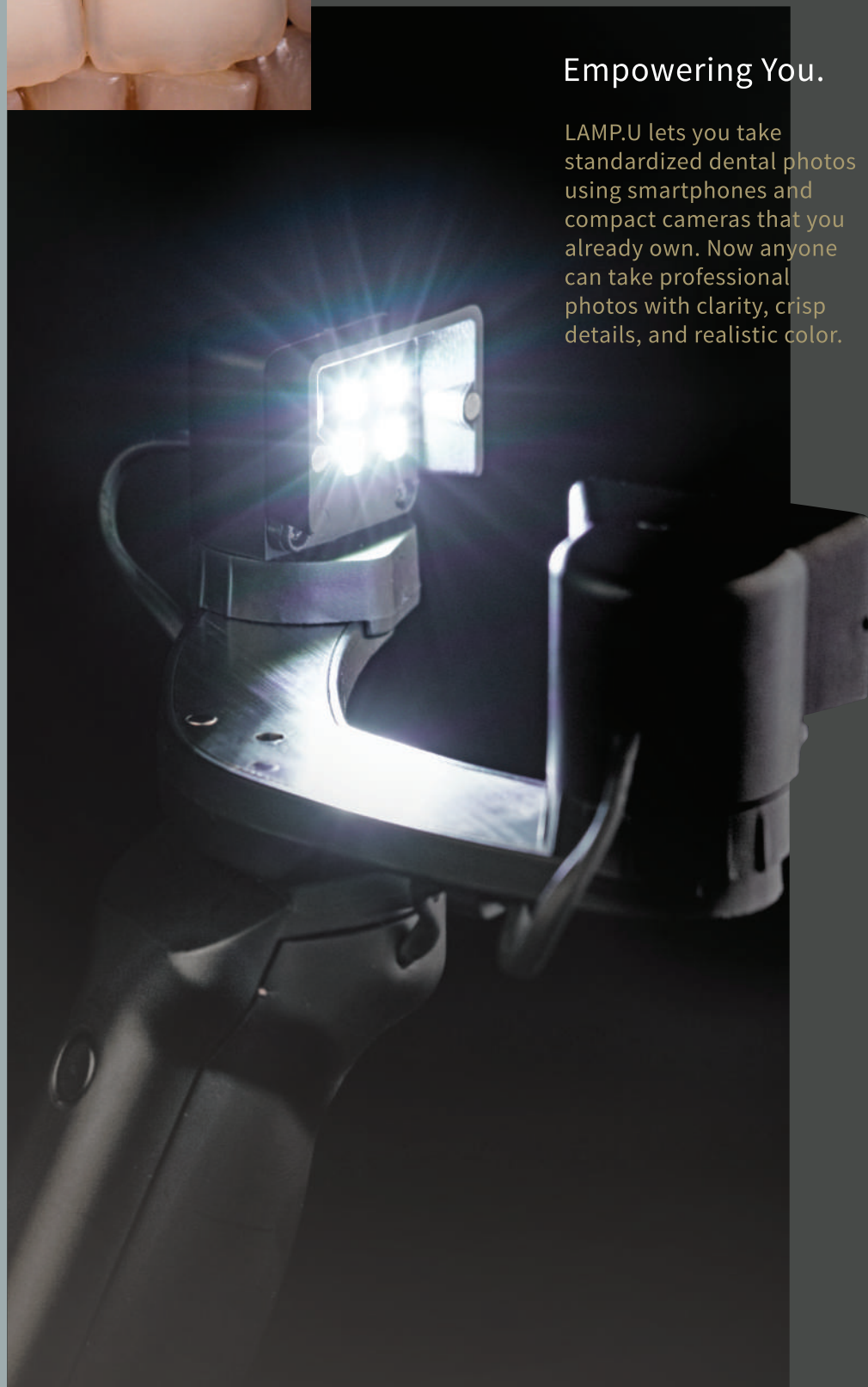




Works Brilliantly on Any Devices.

LAMP.U works really well on a wide range of digital devices — from Android or iOS smartphones to compact digital cameras.

Perfect Images Everytime.

LAMP.U is designed with the beginner dental photographer in mind. Ergonomic to use with its chin guard setup, LAMP.U takes perfect standardized picture - everytime.

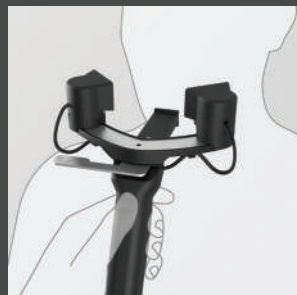

Lateral & Ring LED Light

LAMP.U features two high CRI LED lights which are able to easily slide to two optimized positions. The outer position will allow you to create the perfect 45° angle to light the patient's teeth laterally. The inner position will position the light akin to a ring light, enabling high quality posterior teeth photographs.



95 CRI

75 CRI

High CRI

Bring Out The True Color

CRI (Color Rendering Index) is a measure of how accurately a light source illuminates an object's true color. LAMP.U's LED lights have a CRI rate of up to 98, indicating they are able to produce white light approximating halogen or incandescent lighting and natural daylight.



Cross Polarization Filter

Reflection-Free Images

Three filters snapped into place at LAMP.U's LED and the camera's lens enable teeth photographs without reflection, allowing dentists to accurately check for their true colors.

FEATURES

Lateral and Posterior Light Modes

3 Shooting Modes

- Normal
- Diffused
- Polarized


iOS & Android Compatibility

Digital Camera Compatibility

One Hand System Ergonomic Design

Built-in One Shutter Button via Bluetooth

USB Charging Port

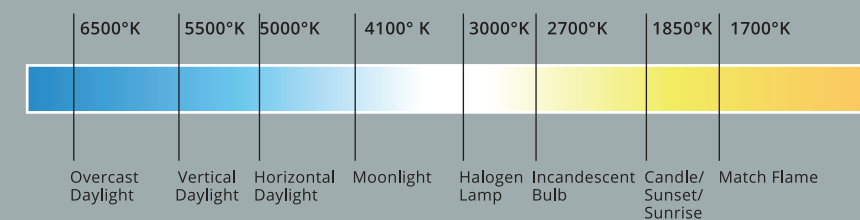

Natural Light

Our LED produces high quality 5500° Kelvin color temperature to imitate the daylight sun, which brings out the teeth's true color.



LAMP.U

Innovative,
User-Friendly and
Affordable.

With LAMP.U, we aim to
help dentists produce
standardized high-quality
dental photographs in
their works.



Distributed by:
Titan Dental Solutions Sdn. Bhd. (1211585-W)

No. 20-2, Jalan 27/70A,
Desa Sri Hartamas,
50480 Kuala Lumpur.
Tel: 03 6206 1175 Sales: 012 731 1529



www.titandental.com.my



[TitanDentalSolutions](https://www.facebook.com/TitanDentalSolutions)

D E N T I T Y



Manufactured by DENTISTRY
in Indonesia

See tutorial at
www.dentistry.net

Contact
info@titandental.com.my

 [@dentistryindonesia](https://www.instagram.com/dentistryindonesia)